




odel x-ray

# INSTALLATION NOTE

## CONNECTION INSTRUCTION NAVIGATOR ON MPU PCB 5806 V4 R8 (and more)


IMI-0167

V1 R1

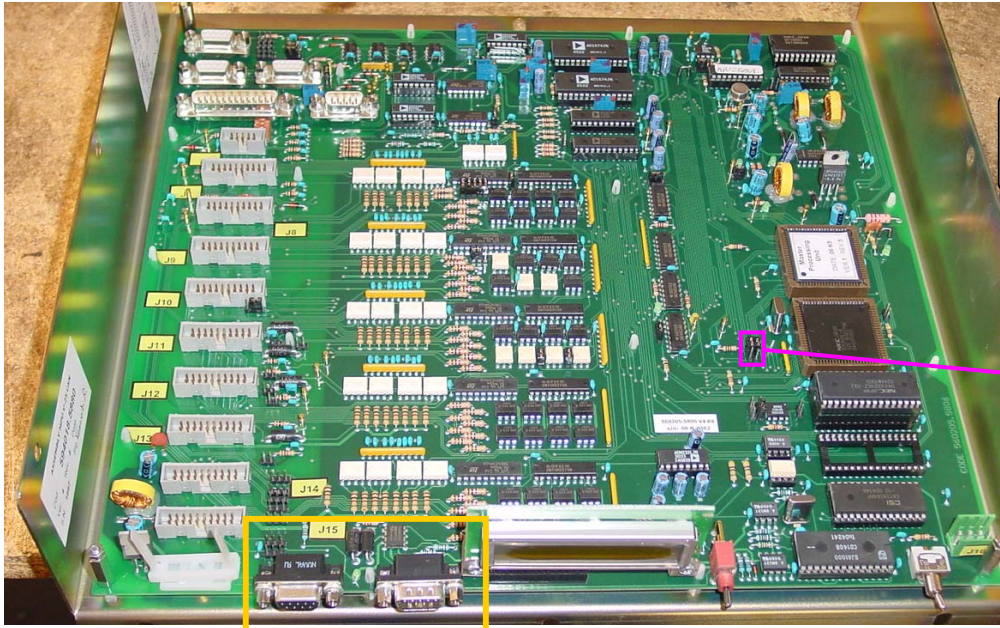

**MCU PCB 5806  
Cod. 560205**

	CONSOLLE	CONNECTION
	NAVIGATOR WITH TOUCH CONSOLLE	V25 CONTROL CONSOLLE
JP100		
JP101		

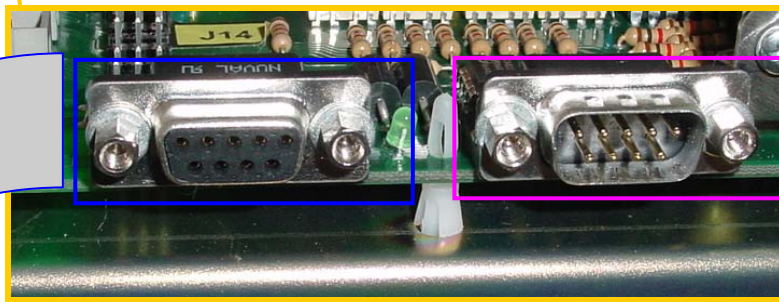

**USE TO CONNECT  
NAVIGATOR DB9  
WITH TOUCH CONSOLLE  
DIRECT CONNECTION**

**USE TO CONNECT  
NAVIGATOR DB9 WITH  
V25 CONTROL CONSOLLE  
NULL MODEM ADAPTOR**

Emesso da: M. Besana	Data Emissione 28/02/06
Approvato da: G. Ganzetti	CLASSE 13
NOME FILE E LOCAZIONE	FOGLIO 1 di 1
ISO 9000/ 13 IMI / IMI-0167 NAVIGATION CONNECTION WITH MPU 5806 V4R8	