



odel x-ray

# INSTALLATION NOTE


## HIRIS RAD TO R306.3X SERIAL LINK

IMI-0170


V1 R1

### *HIRIS RAD TO R306.3X SERIAL LINK*


#### *WITH MALE/FEMALE SERIAL CABLE*

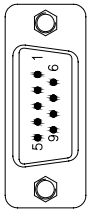


DB9M  
MCU J17

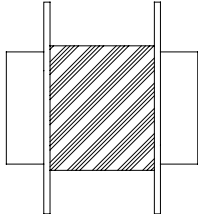


DB9 F/F  
NULL MODEM  
ADAPTER

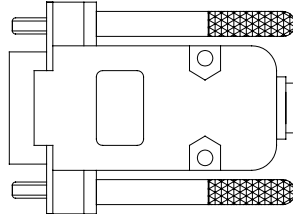

DB9M CABLE  
TO MCU  
RS232 OUT

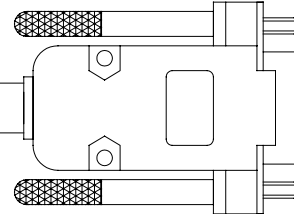


DB9F CABLE  
TO HIRIS  
RS232 OUT

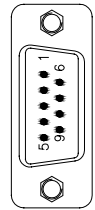


DB9M  
CM24  
HIRIS RAD


#### *WITH FEMALE/FEMALE SERIAL CABLE*

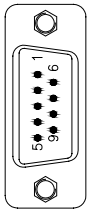


DB9M  
MCU J17

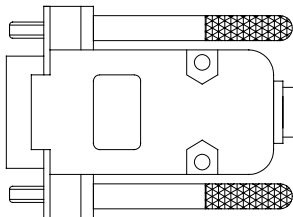

DB9F CABLE  
TO MCU  
RS232 OUT

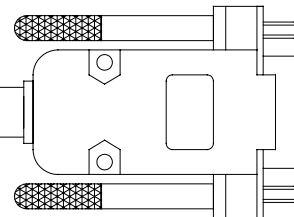

DB9F CABLE  
TO HIRIS  
RS232 OUT

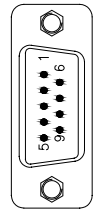

DB9M  
CM24  
HIRIS RAD

Emesso da: <b>G. Bramati</b>	Data Emissione 07/04/2006
Approvato da: <b>S. Ganzetti</b>	CLASSE 13
NOME FILE E LOCAZIONE	FOGLIO 1 di 1
ISO 9000/ 13 IMI / IMI-0170 V1R1 HIRIS RAD TO R3063X SERIAL LINK	